



Clopidogrel is an antiplatelet agent. It is used for treatment of thrombosis, such as stroke and myocardial infarction. Aspirin, which is famous as a non-steroidal anti-inflammatory drug (NSAID), also has an antiplatelet effect.

Advantages and safety of combination drug of clopidogrel and aspirin have been studied and discussed. In this note, clopidogrel and aspirin were successfully separated by using Inertsil ODS-4 with superior peak shapes.

Chemical Structure

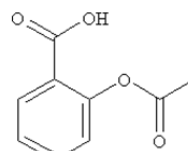


Aspirin

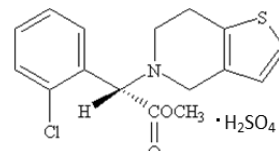

Clopidogrel Sulfate

Structures are created using Chemistry 4-D Draw which is provided by ChemInnovation Software, Inc.

Standard Solution

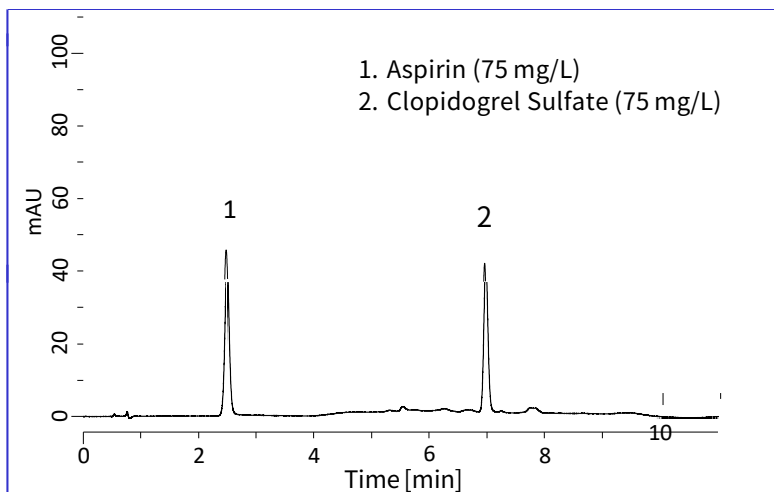


Conditions


Column	: Inertsil ODS-4 HP (3 μm, 50 x 4.6 mm I.D.)
Column Cat. No.	: 5020-14068
Eluent	: A) CH ₃ CN B) 20 mM KH ₂ PO ₄ (pH 3.0, H ₃ PO ₄) A/B = 25/75 -2.5 min- 25/75 -1.5 min- 70/30 -4 min- 70/30 -1 min-25/75 -3 min- 25/75, v/v : 1.0 mL/min
Flow rate	
Column Temp.	: 25 °C
Detection	: UV 240nm
Injection Vol.	: 5 μL

Tablet Sample

The sample concentration was adjusted to 150 mg/L and 75 mg/L as aspirin (acetylsalicylic acid) and clopidogrel, respectively.

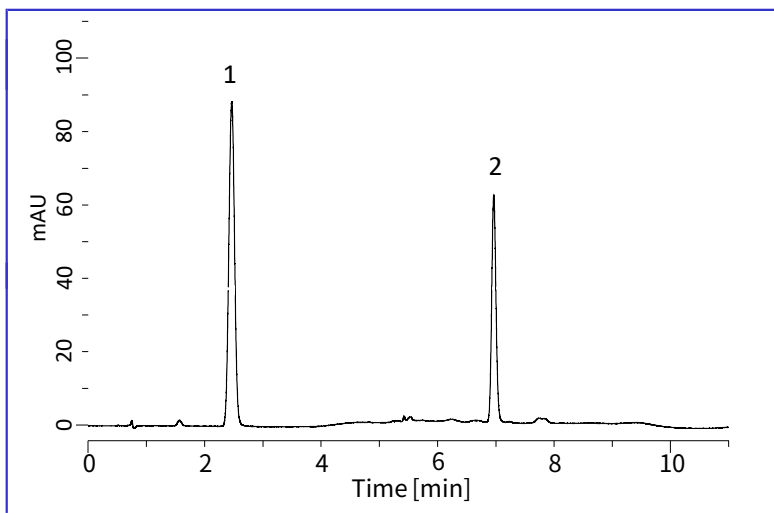

Example of Pretreatment Method for Tablet Sample

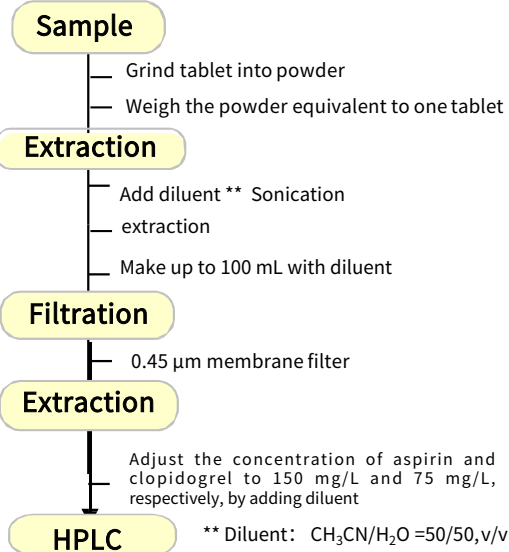

GL Sciences disclaims any and all responsibility for any injury or damage which may be caused by this data directly or indirectly. We reserve the right to amend this information or data at any time and without any prior announcement.

GL Sciences Inc. Japan

22-1 Nishishinjuku 6-chome
Shinjuku-ku, Tokyo
163-1130, Japan

Phone: +81-3-5323-6620
Fax: +81-3-5323-6621
Email: world@glsciences.com
Web: www.glsciences.com

GL Sciences Inc. USA

4733 Torrance Blvd. Suite 255
Torrance, CA 90503
USA

Phone: +1-310-265-4424
Fax: +1-310-265-4425
Email: info@glsciencesinc.com
Web: www.glsciencesinc.com

GL Sciences B.V.

Dillenburgstraat 7C
5652AM, Eindhoven
The Netherlands

Phone: +31-40-254-9531
Email: info@glsciences.eu
Web: www.glsciences.eu

GL Sciences (Shanghai) Limited

Tower B, Room 2003
Far East International Plaza
No.317 Xianxia Road, Changning District
Shanghai, China 200051

Phone: +86-21-62782272
Email: contact@glsciences.com.cn
Web: www.glsciences.com.cn



International Distributors

Visit our Website at www.glsciences.com/distributors